



BOTTLE TOP DISPENSER

Redefining Dispensing
A New era of Precision & Efficiency

Your choice for safe
and precise

Laboratory Dispensing

DIY
DO IT YOURSELF

Re-Calibration

Hybrid valve for smooth operation

Low resistance to
piston movement for
effortless dispensing
of liquid



Easy calibration reset

No need for expensive
third party calibration
after dismantling or
autoclaving

Users can easily reset
the calibration using
the built in tool and
position marker

Easy volume selection

Volume selection is
done by the easy to
read large format
scale with slide &
lock knob

Splash SAFETY Guard

Prevents accidental
splash possibility in
the event of high
pressure down
stroke

Adapter set

Comes with a set of
adapters to mount
on varied bottle
necks




Backed by NABL Lab accreditation

our Bottle Top
Dispenser (BTD)
certificates provide
complete accuracy
and precision
readings, delivering
a new standard of
transparency in the
industry.



Wipeout
piston for
**Smooth
Operation**

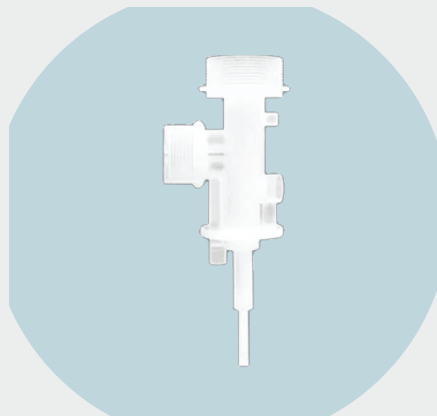

**Springless
Valve**

To ensure a durable,
maintenance free
design



**Safety
Mechanism**

to avoid accidental
push outs of
the discharge valve



**DIY
Calibration**



QR Code
for easy access
to Chemical
Compatibility Chart

Ordering Information



**Bottle Top Dispenser Research Model- Fully Autoclavable,
Calibration and Accuracy as per ISO 8655**

Vol. Range (ml)	Increment (ml)	Accuracy		CV		Product Code RESEARCH Without Recirculation
		± %	± ml	± %	± ml	
0.25 - 2.5	0.05	0.6	0.015	0.2	0.005	BLFLLH002011002D00
0.5 - 5.0	0.1	0.6	0.03	0.2	0.01	BLFLLH002011005D00
1.0 - 10.0	0.2	0.6	0.06	0.2	0.02	BLFLLH002011016D00
2.5 - 25.0	0.5	0.6	0.15	0.2	0.05	BLFLLH002011020D00
5.0 - 50.0	1	0.6	0.30	0.2	0.10	BLFLLH002011040D00
10.0 - 100.0	2	0.6	0.6	0.2	0.2	BLFLLH002011100D00

Dispensers come with 28mm, 32mm, 36mm, 38mm, 40mm,& 45mm bottle neck adapter.



**Bottle Top Dispenser Pinnacle Model- with Re-circulation Valve,
Calibration and Accuracy as per ISO 8655**

Vol. Range (ml)	Increment (ml)	Accuracy		CV		Product Code PINNACLE With Recirculation Valve
		± %	± ml	± %	± ml	
0.25 - 2.5	0.05	0.6	0.015	0.2	0.005	BLFLLH002012002D00
0.5 - 5.0	0.1	0.6	0.03	0.2	0.01	BLFLLH002012005D00
1.0 - 10.0	0.2	0.6	0.06	0.2	0.02	BLFLLH002012016D00
2.5 - 25.0	0.5	0.6	0.15	0.2	0.05	BLFLLH002012020D00
5.0 - 50.0	1	0.6	0.30	0.2	0.10	BLFLLH002012040D00
10.0 - 100.0	2	0.6	0.6	0.2	0.2	BLFLLH002012040D00

Dispensers come with 28mm, 32mm, 36mm, 38mm, 40mm,& 45mm bottle neck adapter.

Note: CV - Coefficient of Variation is also known as Relative Standard Deviation (RSD)

BOROSIL®
Scientific

1101, Crescenzo, G-Block, Opp. MCA Club,
Bandra Kurla Complex, Bandra (E), Mumbai - 400 051, India
E: bsl@borosil.com W: www.borosilscientific.com

Pune Manufacturing Unit:
BOROSIL SCIENTIFIC LTD.
21-24, Plot No. 8, Indialand Global Industrial Park,
Hinjewadi Phase - I, Pune - 411057

Service Support:

✉ lab.support@borosil.com
☎ 1800-224-551

Visit Website



Connect with us
on LinkedIn

